

Shilog Medical Supply

~ ADULT MANUAL OR ELECTRIC BED ~



Model #110E is shown above. Manual version (110M) is also available

Standard Features

- Orthopedically designed four section sleep deck for maximum patient comfort –can be put into the following common positions:
 - Contour
 - Fowler and Semi-Fowler
 - Frame Trendelenburg and Reverse Trendelenburg Positions
 - Cardiac Chair
 - Nursing or Feeding
 - Vascular
- Sleep Deck Dimensions 35" X 80" ~ Weight Capacity is 450 Lbs.
- Sleep Deck Low Height 15.500"
- Sleep Deck High Position 26" (32" with 6" Mattress) Meets Stretcher Height
- Built In Stabilizer Legs (Lifts Bed off Casters in Low Position - Prevents Bed Movement)
- Staff Control Panel at foot end of bed, allows operation of all bed functions as well as lock-out of patient control if desired
- Auxiliary crank for power failure emergencies
- Spring Link Fabric Sleep Deck (Shown with Upgraded Perforated Steel Pan Deck)
- 4" Double Locking Casters (Shown with Deluxe 5" Steer / Total Lock Casters)

- Six Integrated IV Sockets (2 - ½ “ and 2 - ¾” Sockets at head end, 2 – ½” at seat)
- Removable high impact polypropylene foot and head panels with built in tubular push bars for safe and efficient transport.
- Undercarriage is a uni-body welded base for strength and durability.
- All electronic controls are low voltage and UL listed
- Motor drive units are AC powered with 110/115V or 220/240V option and are UL Listed

Options

- Perforated Steel Pan Deck for enhanced support and improved cleanability.
- 5” Deluxe Steer and Total Lock Caster System, allows for easy maneuverability in transport conditions


Side Rails Choices

- Half Head Rail – mounts to movable head section and swings down for entry and exit of patient. (Meets the proposed FDA Guidelines for entrapment)
- Split Rail system – Half Foot Rail is added to the above head rail configuration – adds additional protection for Patient. (Meets the proposed FDA Guidelines for entrapment)
- Three Quarter Length Rail System – offers continuous coverage, 54” length, (Meets the proposed FDA Guidelines for entrapment)




**Above shows the location of the Staff Control
Allows staff to operate bed and lock out the
patient hand control**

**Six function Hand Control allows patient
ability to operate all bed functions**

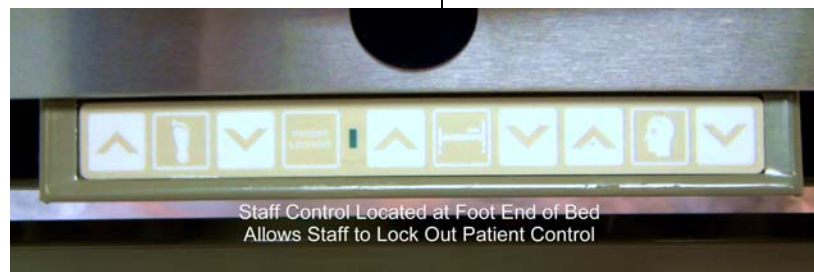


The photo below shows the bed in the Reverse Trendelenburg position

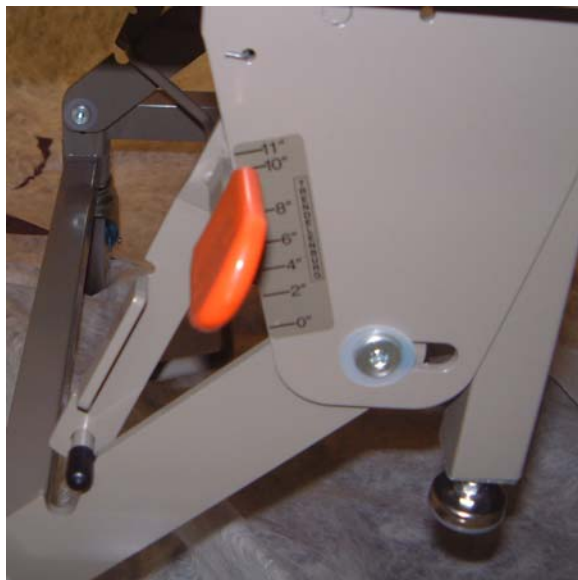

To operate the Trendelenburg and Reverse Trendelenburg, simply raise the bed to the high position and engage the lever as shown. Now press the bed down button to achieve the position.



Stabilizer Legs

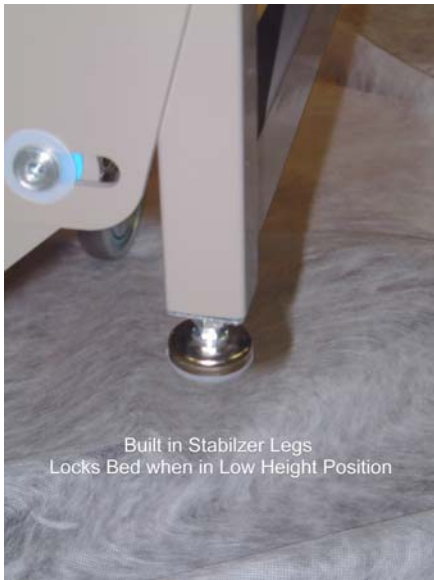

5" Deluxe Caster Option

